




Cleon™

Glide X Operating instructions

Model number: Glide X
Please read this manual carefully before use
The company reserves the right to interpret the specification
The appearance of the product is subject to the actual product
Please keep it together with the invoice after reading
In case of product technology or software upgrade, without prior notice

We focus on providing more scientific and comprehensive oral health care solutions for people. It has been recognized by many consumers. To enable you to use it better, please read this product instructions in detail.

Notes:
Please read all instructions carefully before using the product. If there is any error, please correct us.

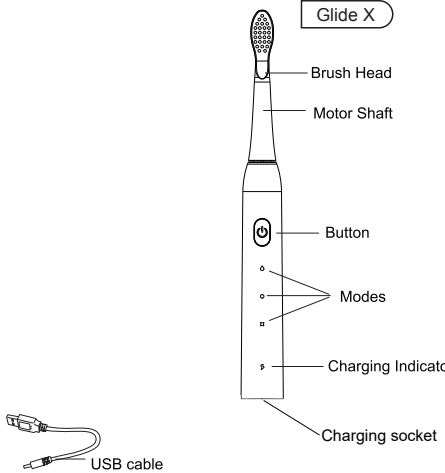

Specification&Information

Item Name	Electric Toothbrush
Model No.	Glide X
Waterproof Level	IPX7
Color	subject to the products you receive

Rated Voltage/Power	USB5.0V === 1.8W
Charging Time	USB charging cable 3 hours, wireless charging 6 hours
Standards	GB4706.1-2005, GB4706.59-2008
Production Batch	Please check the bottom of the product


Packing Accessories

Electric Toothbrush	*1
Brush Heads	*6
Instructions	*1
Charging Cable	*1


Steps:




1 When the power is weak and the LED power indicator flashes red, please charge first. Please refer to the charging instructions for charging.

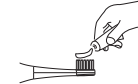

3 Soak the toothbrush head and squeeze an appropriate amount of toothpaste.

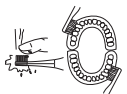

5 Move up and down, do not press hard or scrub.

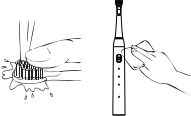


7 After brushing your teeth, please rinse the brush head.




2 Please install the brush head as shown in the picture. There is some gap between the brush head and the handle, which is normal.



4 Put the toothbrush head in your mouth, the bristles and gums are at 45°, activate the switch, select the right brushing mode for you.

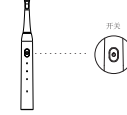

6 To prevent the toothpaste foam from splashing, please take it out of your mouth after turning off the.

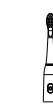

8 Place brush head in ventilated dry place

Brushing Strengthsthe

Glide X offers three different strengths forvarying oral care needs:

Soft - For sensitive teeth and gums, brushingyour teeth in gentle frequency.

Medium - Switching continuously and fast andwhitening teeth.

Strong - Brush head swing in high frequency,remove tartar and bacteria on teeth.

Disposal

Important: do not clean the charging unit under running water or immerse it in water.

Place the charging station on dry and safe place.

This product contains a rechargeable battery. In the interest of protection the environment, please do not dispose of the product in the household waste at the end of its useful life.

Charging Description

- Please check the power supply voltage and the power socket are safe and normal, and the voltage adaption.
- Please make sure that the USB cable is dry, the charging plug is free from water stains and damage, etc., and then plug in the power supply under safe conditions.
- When using an electric toothbrush for the first time, please turn on the electric toothbrush first. In case of low or no power, please fully charge the electric toothbrush first.

Battery Indicator Description

Indicator status during discharge	Description
Indicator on	Standby Modes
Red indicator on	Low battery, please charge in time
Indicator not on	Ultra low battery, products cannot work normally
Indicator status during charging	Description
Red indicator on	on charging
Indicator turns from red to green	Finished charging

Notes

- Please pay attention to the status of the battery indicator during use.
 - Please charge in time when the battery is low
- If you have not used the electric toothbrush for a long time, please charge the toothbrush when you use it again. And to ensure safe charging, which is conducive to improving the service life of the product.



Warantee Card

Item	
Model No.	
Date of Purchase	
Name	
Mobile	
Address	
Fault Description	

Please fill in your information carefully, and please send this card together with the product to be repaired to our repair center.

For Any Feedback / Complaints :

+91 9979175505
care@cleonsmart.com

Dental Expert Advice

Brush your teeth not less than twice a day.
Replace the brush head every 2-3 months.

Cleon™