

Cleon™



Orvix Operating instructions

* product introduction

This electric toothbrush is a high-frequency sweeping electric toothbrush with a vibration frequency of tens of thousands of times. The high-frequency sweeping vibration can easily remove germs and impurities around the teeth, and its cleaning range can cover all parts around the teeth. The acoustic energy generated by the motor can be transmitted to the surfaces of teeth and gums through the bristles of the brush head, so as to remove tartar and bacteria.


* Attention

1. Please read this user manual carefully before using the product. This product is suitable for children and adults. Children should use it under the supervision of adults and follow the following basic safety precautions.
2. Users of orthodontics, dental prostheses and dental implants should use them under the guidance of dentists.
3. Please consult the doctor or the manufacturer of the implanted device before using any questions related to the pacemaker or other implanted devices.
4. This product is only designed for cleaning teeth, gums and tongue. Do not use it for other purposes. If you feel discomfort or pain, please stop using this product and contact your doctor.
5. Please stop using the brush head with deformed or bent bristles. Please replace the toothbrush head every three months. If there are signs of wear, please replace it in advance. Do not use other brush heads that are not recommended by the manufacturer.
6. If your toothpaste contains peroxide, sodium bicarbonate or other bicarbonate (common ingredients in whitening toothpaste), please be sure to thoroughly clean the brush head with fertilizer and water after each use to prevent the plastic from cracking.

* Danger

1. Please check whether the product is obviously damaged before each use. Please stop using the machine if it is damaged or seriously dropped.
2. If the power cord or plug (power cord) is damaged, please be sure to replace it with the original model.
3. If this product is damaged (brush head and handle), please stop using it.
4. It is forbidden to charge in the bathroom or in a humid environment; Please do not charge this product for a long time to avoid unnecessary danger.

Product description

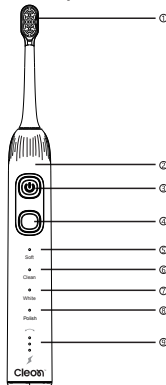

- ① Toothbrush head
- ② Toothbrush handle
- ③ Switch and mode adjustment buttons
- ④ Force keying
- ⑤ Softing mode
- ⑥ Cleaning mode
- ⑦ Whitening mode
- ⑧ Polishing mode
- ⑨ Force indicator light
- ⑩ Charging cable

* Warning

1. Please check whether the product is obviously damaged before each use. Please stop using the machine if it is damaged or seriously dropped.
2. If the power cord or plug (power cord) is damaged, please be sure to replace it with the original model.
3. If this product is damaged (brush head and handle), please stop using it.
4. It is forbidden to charge in the bathroom or in a humid environment; Please do not charge this product for a long time to avoid unnecessary danger.

* Reminder

1. Do not share the brush head with others.
2. Do not wrap the power cord around the adapter when storing.
3. Do not put metal objects, such as coins or paper clips, on the adapter.
4. If it is not suitable for a long time, please unplug the power supply to avoid fire.

* Function introduction

1. This product has a one-button start system, including on/off and gear shift. 4 brushing modes, each mode has 3 strength choices.

2. Before use, press the upper button to start the machine, select the brushing mode, and then select the brushing intensity through the lower button.
3. When each mode is working, pause for 0.3 seconds every 30 seconds, press stop or restart halfway, and the working time will be counted again.

* mode introduction:

- Soft sweep vibration mode--Soft vibration sweeping mode, three kinds of force adjustment, soft gum protection, sensitive for novices. (2 minutes).
- Clean sweep vibration mode--Cleaning vibration sweeping mode, three kinds of force adjustment, deep cleaning between teeth (2 minutes).
- White sweep vibration mode--Strong vibration sweeping mode, three kinds of force adjustment, deep stain removal, clean and worry-free teeth. (2 minutes).
- Polishing mode -Polish mode, three kinds of intensity adjustment, net removal of smoke stains and tartar (2 minutes).

* Charging

1. When working, the power indicator light is green. When the (3.3V) indicator light changes from green to red and flashes, the main body needs to be charged. This product has the function of battery over-discharge protection. When the battery is not charged at low voltage in time, it will be forced to shut down and cannot be used normally.
2. Charging method: Type-C charging cable charging.
3. The toothbrush can be charged for about 8 hours. When charging, just connect the charging cable to the power supply.
4. Charging: When charging, the indicator light is red and flashes. When fully charged, it turns from red to green and stays on.
5. If the product is left for a long time without charging, the indicator light may not light up immediately when charging at this time, which is normal. Generally, it takes 30 minutes to pre-charge the product before the indicator light can work normally.

* Preparation for use

1. Align the front section of the toothbrush head with the front section of the toothbrush handle.
2. The toothbrush head is installed on the metal shaft (there will be a small gap between the two parts, so that the toothbrush head does not collide with the main engine during work, so that the motor has stronger power).

* Brush your teeth

Step 1

Wet the bristles and squeeze some toothpaste on the brush head.

Step 2

Press the on/off button to open the toothbrush.

Move the toothbrush slowly at an angle of 45 between the bristles and the gums.

Press the button during brushing to switch the brushing mode freely.

It is recommended to brush your teeth for 2 minutes.

* clean the brush head

After brushing your teeth, aim the toothbrush head at the water and clean it. After cleaning, cover the brush head to protect the brush head clean and hygienic.

* Maintenance instructions

1. Wash the brush head and bristles after each use and air dry naturally. (Please don't pull the bristles, it will easily deform them and fall off.
 2. Wipe with a dry cloth after each use and air dry naturally.
 3. Do not clean the main body, power cord and plug with water above 40°C or acidic cleaning solution. It is forbidden to immerse the plug in water.
 4. When it is not charged, be sure to pull out the power plug from the socket, otherwise it will easily lead to electric shock.
 5. It is recommended to replace the brush head every three months, or when the bristles are bent or deformed.
- * Environmental protection**
- When disposing of products, please don't pile them together with general domestic garbage.
 - The built-in rechargeable battery contains substances that may pollute the environment. Please don't throw away the product at will. Please put the product in front of the recycling bin designated by the government.

* Specification parameters

Product model: Orvix
Battery voltage: 3.7V
Rated input: 5V 1A
Rated power: 5W
Waterproof grade: IPX7
Battery capacity: LI-ION 1200MAH
Vibration frequency: 30000VPM

For Any Feedback / Complaints :

+91 9979175505
care@cleonsmart.com



1. It is forbidden to disassemble, impact, squeeze or throw into the fire.
2. If there is serious swelling, please do not continue to use it.
3. Do not place in high temperature environment.
4. Keep the charging port dry when charging.
5. It is forbidden to use the battery after it is flooded.
6. Before discarding the appliance, the battery must be removed from the appliance.
7. When taking out the battery, the appliance must be powered off.
8. The battery should be disposed of safely.

* Description of hazardous substances

Name and content of harmful substances in products

Component name	Harmful substance					
	lead (Pb)	mercury (Hg)	cadmium (Cd)	hexavalent chrome (Cr(VI))	polybrominated biphenyl (PBB)	polybrominated diphenyl ether (PBDE)
Brush head	○	○	○	○	○	○
Toothbrush handle	○	○	○	○	○	○
Electric cell	○	○	○	○	○	○
Attachment	○	○	○	○	○	○

This form is compiled according to SJ/T11364.

- : It is indicated that the content of the harmful substance in all homogeneous materials of this part is within the limit specified in GB/T26572.
- x : indicates that the content of the harmful substance in at least one homogeneous material of the component exceeds the limit requirements specified in GB/T26572.

* Warranty card

Thank you for your support and trust!

The Orvix toothbrush should be able to give full play to its advantages in normal use and maintenance. Any damage caused by manufacturing process or components within one year of purchase will be replaced free of charge.

* Warranty terms

1. Free warranty date, starting from the date of purchase.
2. Please show the back voucher or purchase record for the warranty products.
3. Replacement of spare parts for more than one year will be charged according to the replacement price.
4. The toothbrush head is a consumable product, which is not covered by the warranty.
5. Under any of the following circumstances, free warranty is not implemented.
 - 1) Dismantle the body by yourself during use, resulting in the toothbrush not being used normally.
 - 2) Unable to provide purchase voucher or purchase record.