


Cleon™

Swypon Operating instructions

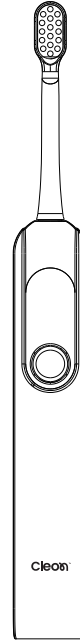

* product introduction

This electric toothbrush is a high-frequency sweeping electric toothbrush with a vibration frequency of tens of thousands of times. The high-frequency sweeping vibration can easily remove germs and impurities around the teeth, and its cleaning range can cover all parts around the teeth. The acoustic energy generated by the motor can be transmitted to the surfaces of teeth and gums through the bristles of the brush head, so as to remove tartar and bacteria.


* Attention

1. Please read this user manual carefully before using the product. This product is suitable for children and adults. Children should use it under the supervision of adults and follow the following basic safety precautions.
2. Users of orthodontics, dental prostheses and dental implants should use them under the guidance of dentists.
3. Please consult the doctor or the manufacturer of the implanted device before using any questions related to the pacemaker or other implanted devices.
4. This product is only designed for cleaning teeth, gums and tongue. Do not use it for other purposes. If you feel discomfort or pain, please stop using this product and contact your doctor.
5. Please stop using the brush head with deformed or bent bristles. Please replace the toothbrush head every three months. If there are signs of wear, please replace it in advance. Do not use other brush heads that are not recommended by the manufacturer.
6. If your toothpaste contains peroxide, sodium bicarbonate or other bicarbonate (common ingredients in whitening toothpaste), please be sure to thoroughly clean the brush head with fertilizer and water after each use to prevent the plastic from cracking.

* Danger

1. Please check whether the product is obviously damaged before each use. Please stop using the machine if it is damaged or seriously dropped.
2. If the power cord or plug (power cord) is damaged, please be sure to replace it with the original model.
3. If this product is damaged (brush head and handle), please stop using it.
4. It is forbidden to charge in the bathroom or in a humid environment; Please do not charge this product for a long time to avoid unnecessary danger.

Product description

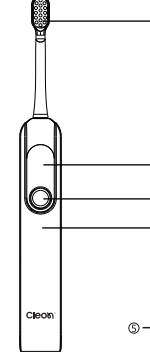

- ① Toothbrush head
- ② Toothbrush host
- ③ Mode and function display
- ④ Switch and mode adjustment buttons
- ⑤ Charging cable

* Warning

1. Please check whether the product is obviously damaged before each use. Please stop using the machine if it is damaged or seriously dropped.
2. If the power cord or plug (power cord) is damaged, please be sure to replace it with the original model.
3. If this product is damaged (brush head and handle), please stop using it.
4. It is forbidden to charge in the bathroom or in a humid environment; Please do not charge this product for a long time to avoid unnecessary danger.

* Reminder

1. Do not share the brush head with others.
2. Do not wrap the power cord around the adapter when storing.
3. Do not put metal objects, such as coins or paper clips, on the adapter.
4. If it is not suitable for a long time, please unplug the power supply to avoid fire.

* Function introduction

1. This product has a one-button start system, including on/off and gear shift.

2. Press the power-on button briefly, enter the mode interface after 3 seconds, press the button to switch the mode within 5 seconds, and press the button to turn it off after 5 seconds.
3. When each mode is working, pause for 0.3 seconds every 30 seconds, press stop or restart halfway, and the working time will be counted again.
4. In the power-off state, press and hold the power-on button for more than 5 seconds to enter the reset system interface, and then press it again to select languages. In the selected language interface, press and hold the power-on button for more than 2 seconds to select the language successfully.

* mode introduction:

- Soft: A soft vibration frequency designed for first-time users of sonic toothbrushes (2 minutes);
- Gentle: A comfortable vibration frequency for sensitive teeth and gums (2 minutes);
- Clean: High-speed vibration frequency bristles can easily remove tartar and destroy bacteria on the tooth surface and in hidden areas (2 minutes);
- White: An ultra-high frequency sonic mode that can polish and brighten teeth (2 minutes).

* Charging

1. When charging, the battery logo on the screen shows that it is charging, and when it is fully charged, it shows that it is full.
2. When the power is less than 10%, please charge it in time; when it is less than 1%, it cannot be turned on.
3. charging method: Type-C charging line charging.
4. The toothbrush can be charged for about 7 hours. When charging, just connect the charging cable to the power supply.
5. If the product is left for a long time without charging, the screen may not light up immediately when charging, which is normal. Generally, it takes 30 minutes to pre-charge before the screen can work normally.

* Preparation for use

1. Align the front section of the toothbrush head with the front section of the toothbrush handle.
2. The toothbrush head is installed on the metal shaft (there will be a small gap between the two parts, so that the toothbrush head does not collide with the main engine during work, so that the motor has stronger power).

Function Introduction:

Power-on function	1. In the awakened state, a short press will trigger and enter the power-on and gear-switching state, while the battery level will be displayed. 2. For the first power-on, the default gear position is set to one. For subsequent power-ons, the previous working gear position will be remembered.
Language Settings	In the awakened state, press the button for 3 consecutive short times to enter the language design menu. The details can be found in "Language Settings". Press the switch lightly to select the language design, and press and hold for 2 seconds to exit and select the current language. The default language is English. The language options are English, Chinese, and Korean.
Gearshift	1. When the product is powered on and enters the gear shifting state, press the button within 8 seconds to switch gears. The gear shifting sequence is: gentle → sensitive → clean → whitening → gentle (in a cycle). After shifting gears, the 8-second timer is reset and starts counting again. After 8 seconds, press the power-off button briefly. 2. After stopping the gear shifting for 8 seconds, the operation menu will be entered.
WorkMenu	1. The entire operation time is 2 minutes. A 0.5-second pause prompt is given every 30 seconds. 2. The display time of the content and gear information is provided.
Pick up the wake-up call.	1. When the travel lock is turned off, the toothbrush enters the wake-up state and displays the previous brushing area prompt interface for 5 seconds. 2. If no button is pressed within 5 seconds, it goes into sleep mode; if a short press is made, it will turn on.
Pressure sensing	1. When the force is less than 40gf, the toothbrush enters the anti-splash mode and the motor vibrates weakly. 2. When the force is between 40gf and 300gf, the motor vibrates at the preset speed. 3. When the force is greater than 300gf, an abnormal vibration alert is triggered (vibration switch) and the screen displays an overpressure indicator. The pressure is reduced and the previous working menu is immediately returned. 4. The anti-splash mode is not counted in the brushing time. If the anti-splash mode lasts for more than 30 seconds, it enters the end state.
Travel lock	1. When the device is powered off or in the "not powered on" state, press a and hold for more than 2 seconds to activate the travel lock. At this point, the wake-up function will be deactivated; 2. Press briefly (< 2 seconds) and it will prompt "Press again to deactivate the travel lock". No other actions are required. Press for more than 2 seconds to deactivate the travel lock and then pick up the wake-up function to activate it.

Final state	1. Once the working mode is completed, either by long-pressing for 1 second to force a shutdown, or when the anti-dust state lasts for more than 30 seconds, the system will directly enter the end state. 2. First, an interface for timing and scoring will appear: it will display the coverage of the brushing area and the total brushing duration. After 3 seconds, the interface for area re-brushing will appear. 3. Area prompt interface: The entire mouth is divided into 6 areas: upper left, lower left, upper middle, lower middle, upper right, and lower right. The effective brushing duration in each area exceeding 10 seconds is considered a full score. The coverage score is calculated based on the effective brushing duration in each area. The re-brushing interface shows 10 seconds, and after 10 seconds, it enters the sleep mode. A short press within 10 seconds can restart the device and continue brushing.
Battery level indicator	1. The device will not power on if the voltage is below 3.0V. Each time the power-on button is pressed, a low voltage symbol (lasting for 2 seconds) will be displayed, and it will not start. When the voltage is below 3.0V, the wake-up function will be deactivated until the charging action is detected, then it will be activated. 2. During use: If the voltage is above 3.4V, the icon shows one bar of battery charge, indicating 20% of battery capacity; if the voltage is above 3.6V, the icon shows two bars of battery charge, indicating 40% of battery capacity; if the voltage is above 3.75V, the icon shows three bars of battery charge, indicating 60% of battery capacity; if the voltage is above 3.85V, the icon shows four bars of battery charge, indicating 80% of battery capacity; if the voltage is 4.0V or above, the icon shows five bars of battery charge, indicating 100% of battery capacity. 3. Charging process: If the voltage is less than 3.4V, the icon shows one bar of battery charge (flashing); if the voltage is above 3.4V, the icon shows two bars of battery charge (flashing), indicating 20% of battery capacity; if the voltage is above 3.6V, the icon shows three bars of battery charge (flashing), indicating 40% of battery capacity; if the voltage is above 3.75V, the icon shows four bars of battery charge (flashing), indicating 60% of battery capacity; if the voltage is above 3.85V, the icon shows five bars of battery charge (flashing), indicating 80% of battery capacity; if the voltage is 4.0V or above, the icon shows five bars of battery charge, not flashing, indicating 100% of battery capacity.

* Brush your teeth

Step 1

Wet the bristles and squeeze some toothpaste on the brush head.

Step 2

Press the on/off button to open the toothbrush.

Move the toothbrush slowly at an angle of 45 between the bristles and the gums.

Press the button during brushing to switch the brushing mode freely.

It is recommended to brush your teeth for 2 minutes.

* clean the brush head

After brushing your teeth, aim the toothbrush head at the water and clean it. After cleaning, cover the brush head to protect the brush head clean and hygienic.

* Maintenance instructions

1. Wash the brush head and bristles after each use and air dry naturally. (Please don't pull the bristles, it will easily deform them and fall off.
 2. Wipe with a dry cloth after each use and air dry naturally.
 3. Do not clean the main body, power cord and plug with water above 40°C or acidic cleaning solution. It is forbidden to immerse the plug in water.
 4. When it is not charged, be sure to pull out the power plug from the socket, otherwise it will easily lead to electric shock.
 5. It is recommended to replace the brush head every three months, or when the bristles are bent or deformed.
- * Environmental protection
- When disposing of products, please don't pile them together with general domestic garbage.
- The built-in rechargeable battery contains substances that may pollute the environment. Please don't throw away the product at will. Please put the product in front of the recycling bin designated by the government.

* Specification parameters

Product model: Swypon

Battery: 3.7V

Input voltage: 5V=1A

Power: 2W

Waterproof grade: IPX7

Battery capacity: 1200mAh

Vibration frequency: 38400 Vibrations Per Minute

For Any Feedback / Complaints :

+91 9979175505
care@cleonsmart.com

- ⚠

 1. It is forbidden to disassemble, impact, squeeze or throw into the fire.
 2. If there is serious swelling, please do not continue to use it.
 3. Do not place in high temperature environment.
 4. Keep the charging port dry when charging.
 5. It is forbidden to use the battery after it is flooded.
 6. Before discarding the appliance, the battery must be removed from the appliance.
 7. When taking out the battery, the appliance must be powered off.
 8. The battery should be disposed of safely.

* Description of hazardous substances

Name and content of harmful substances in products

Component name	Harmful substance					
	lead (Pb)	mercury (Hg)	cadmium (Cd)	hexavalent chromium (Cr(VI))	polybrominated biphenyl (PBB)	polybrominated diphenyl ether (PBDE)
Brush head	○	○	○	○	○	○
Toothbrush handle	○	○	○	○	○	○
Electric cell	○	○	○	○	○	○
Attachment	○	○	○	○	○	○

This form is compiled according to SJ/T11364.

○: It is indicated that the content of the harmful substance in all homogeneous materials of this part is within the limit specified in GB/T26572.

×: indicates that the content of the harmful substance in at least one homogeneous material of the component exceeds the limit requirements specified in GB/T26572.

* Warranty card

Thank you for your support and trust!

The electric toothbrush should be able to give full play to its advantages in normal use and maintenance. Any damage caused by manufacturing process or components within one year of purchase will be replaced free of charge.

* Warranty terms

1. Free warranty date, starting from the date of purchase.
2. Please show the back voucher or purchase record for the warranty products.
3. Replacement of spare parts for more than one year will be charged according to the replacement price.
4. The toothbrush head is a consumable product, which is not covered by the warranty.
5. Under any of the following circumstances, free warranty is not implemented.
 - 1) Dismantle the body by yourself during use, resulting in the toothbrush not being used normally.
 - 2) Unable to provide purchase voucher or purchase record.